


# WINDTUNNEL AUTOSAMPLER ENCLOSURE

Protects samples against airborne contaminants  
Protects the operator against fumes and odours

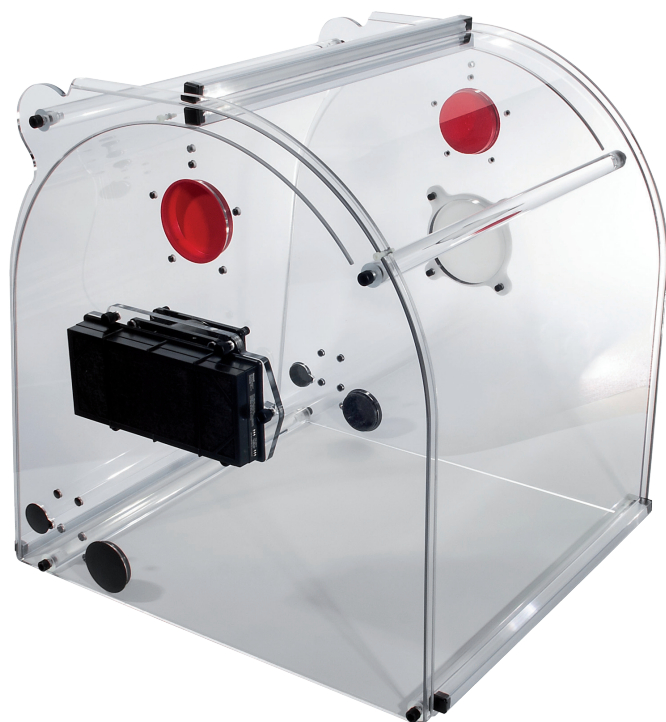

- Suitable for the Cetac 130, 260, 520 and 520HS autosamplers
- Shatter-proof polycarbonate construction
- Easy access through front and rear roll-up doors
- Optional fan maintains positive pressure
- Optional HEPA filter
- Optional exhaust port flange
- Multiple ports for tubing and power cords

## Versatile Configuration

The optional accessories allow you to configure the WindTunnel to best suit your requirements.

**Filtered positive pressure enclosure.** When fitted with the optional fan and HEPA filter, the WindTunnel provides a positive pressure enclosure.

**Vented enclosure.** The addition of the optional flange, allows the fumes to be removed through a laboratory exhaust system (not supplied).

**Dust cover.** The WindTunnel can be used without the fan and filter to function as a convenient dust cover.

The WindTunnel is designed to fit over an autosampler sitting on a solid flat surface. If required, a solid polyethylene base is available separately.

## INTERNATIONAL

Glass Expansion • 6 Central Boulevard • Port Melbourne • Vic 3207, Australia  
Telephone: +61 3 9320 1111 • Toll Free (Aust): 1800 777 638  
Facsimile: +61 3 9320 1112 • Email: [enquiries@geicp.com](mailto:enquiries@geicp.com)

## AMERICAS

Glass Expansion • 4 Barlows Landing Road • Unit 2A • Pocasset, MA 02559, USA  
Telephone: 508 563 1800 • Toll Free (US): 800 208 0097  
Facsimile: 508 563 1802 • Email: [geusa@geicp.com](mailto:geusa@geicp.com)

[www.geicp.com](http://www.geicp.com)



**GLASS EXPANSION**  
Quality By Design